


New Objective Digital PicoView® (DPV) 450 Nanospray Source For Sciex TripleToF 5600 series mass spectrometer

<http://www.newobjective.com/products/sources/dpv-4.shtml>

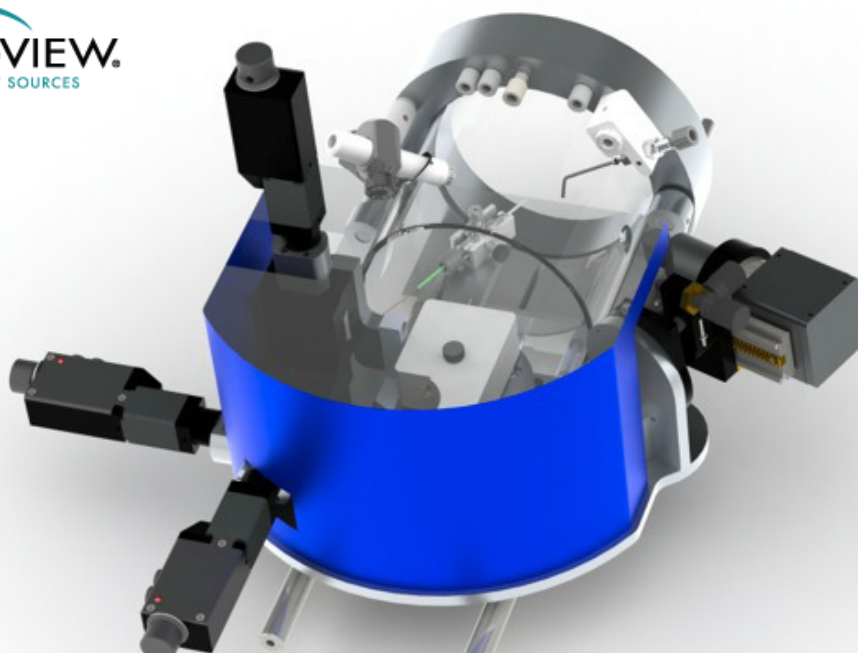

Digital PicoView 450 – What it does for you

- Digitally-controlled fine-precision XYZ stage (1µm precision)

Control, also over remote desktop through the internet

- Ultra-high resolution digital camera

Excellent video spray image, also over remote desktop through the internet

- PV Acquire™ stage control software

Full motorized control over spray positioning

- Automated tip rinsing station

Allows washing of tip between runs and during idle time, thus less clogging and better uptime

- Tool-free magnetic interlock stage plate

Quick swap between different plumbing setups

- Column heating option (single column only)

Less backpressure over room temperature operation – enables longer columns with improved peak capacity

- Dual column option

May double your through put, load and wash channel A while you run channel B

- Sheath gas module included

Counter the curtain gas difficult to handle in NanoESI – larger sweet spot

- Fully compatible with all [PicoTip](#)[®] emitters, [IntegraFrit](#)[™] and [PicoFrit](#)[®] columns

www.mswil.com